



OPT 80 OPU TROLLEY



- ⦿ Designed to support OPU procedures, the OPT 80 enables easy transfer of oocytes into 37°C pre-warmed test tubes by positioning them close to the gynecological table - reducing oocyte stress.
- ⦿ Supports the OPU procedure with its drawer system, which can move up and down and rotate 360° around its axis, offers flexibility and ease of access.
- ⦿ Thanks to its generously sized shelves, the OPT 80 allows the aspiration pump to be used from either side. The three drawers provide convenient storage for consumables both during and after the procedure.
- ⦿ The heated surface, open shelf, and movable mechanism are made from antibacterial Corian material.



CUSTOMIZABLE USE

The OPT 80's movable mechanism tower is designed to be positioned on either the right or left side, depending on your workspace layout and usage preferences.



HEATED SURFACE AREA


With a 450mm x 350mm heated surface, the OPT 80 offers a spacious working area while ensuring stable and efficient heat distribution thanks to its Corian material.






EASE OF USE

Its extendable mechanical design ensures an ergonomic working area, making OPU procedures easier for the team.

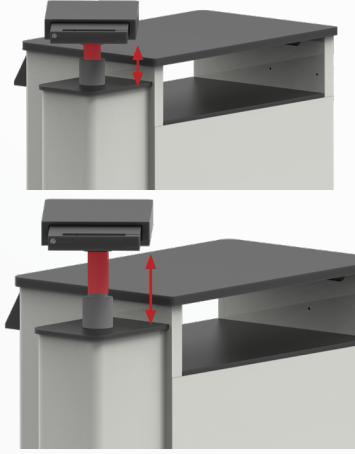

360° SWIVEL HEAD

The mechanical head, which moves up and down and rotates 360°, offers full accessibility from every angle - making it an essential tool in OPU operating rooms for comfortable and efficient procedures

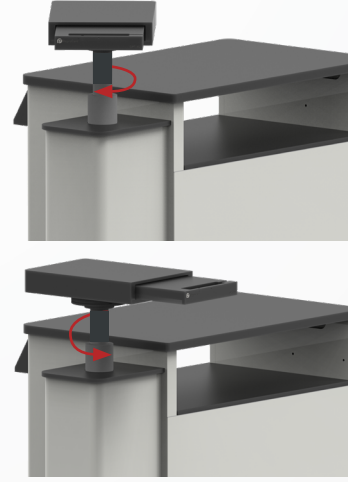

TECHNICAL SPECIFICATIONS

Power Input	220V 50-60Hz
External Dimensions (W x D x H)	1065 x 560 x 820mm
Heated Surface Dimensions (W x D)	450 x 350mm
Weight (incl. Power Adapter)	80kg
Temperature Accuracy	±0.3 °C
Temperature Resolution	0.1 °C
Maximum Heating Power	75 Watt
Nominal Heating Power	55 Watt
Control System	PID Control
Display	4.3" Color Touch TFT
Sensor Type	PT100 Class A

Tekservis Teknik ve Tıbbi Ürünleri Pazarlama San. ve Tic. Ltd. Şti.

Head Office:

Huseyin Rahmi Sk. No: 7 / 7
Cankaya / Ankara / TURKEY

T: +90 (312) 440 7091
F: +90 (312) 440 7090

Branch:

Goztepe Mahallesi,
Goksu Evleri, Hurrem Sokak
No: B248A Beykoz / Istanbul / TURKEY

T: +90 (216) 474 3664
F: +90 (216) 474 0171

www.reprocare.com.tr
info@reprocare.com.tr